

MARQUARDT MEDICAL FACE MASK



KEY FEATURES

- CE certification
- EN14683:2019+AC:2019 – TYPE IIR medical face mask
- Integrated, adjustable and durable nose clip
- Straps tested for skin compatibility

MATERIAL COMPOSITION

Outer and inner non-woven fabrics:

- PP spun-bonded fabric
- PET / PP non-woven fabric for increased wearing comfort

Filter layer:

- PP melt-blown non-woven fabric

Nose clip:

- HDPE-coated galvanized wire, wire thickness: 0.5 mm

Elastic strap:

- Polyamide / elastane elastic straps

INTENDED USE

The medical face mask is designed to protect the patient from foreign bodies, by covering the mouth and nose and providing a barrier to minimize the direct transmission of infectious droplets between staff and patients.

WARNINGS AND LIMITATIONS


- Ensure that you put on the product correctly.
- Dispose the mask at the end of a shift. Do not touch the mask with your hands.
- Wash or disinfect your hands after taking off the mask.
- Wearing the mask does not make you immune to infectious germs.
- Do not use the mask in a way that is not specified in the instructions.

APPROVALS AND STANDARDS

- Medical filtering half-mask according to EN 14683, Type IIR
- Bacterial filtration efficiency $\geq 98\%$



Product number 7112.5001

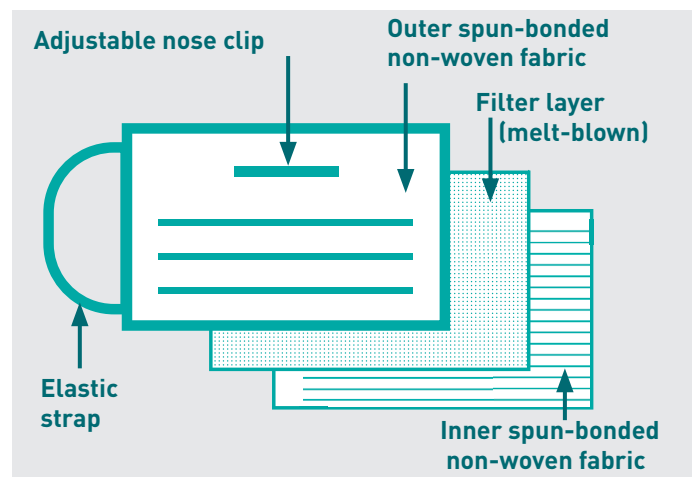

TIME LIMIT OF USE

Single use, max. one shift. Respirator may be used until damages, breathing becomes difficult or the material is soaked with liquid.

SHELF LIFE AND STORAGE

- 2 years from the date of manufacturing in the original packaging.
- Ensure that the product is still within its storage period before use. The product should be stored in a clean and dry place within the temperature range of $-5\text{ }^{\circ}\text{C}$ and $+40\text{ }^{\circ}\text{C}$ at a maximum relative humidity of 60 %.
- Use the original packaging for storage and transport

For further information please refer to the instructions at:

www.marquardt.com/us/masks/

