

# MARQUARDT MEDICAL FACE MASK



## KEY FEATURES

- CE certification
- EN14683:2019+AC:2019 – TYPE IIR medical face mask
- Integrated, adjustable and durable nose clip
- Straps tested for skin compatibility

## MATERIAL COMPOSITION

### Outer and inner non-woven fabrics:

- PP spun-bonded fabric
- PET / PP non-woven fabric for increased wearing comfort

### Filter layer:

- PP melt-blown non-woven fabric

### Nose clip:

- HDPE-coated galvanized wire, wire thickness: 0.5 mm

### Elastic strap:

- Polyamide / elastane elastic straps

## INTENDED USE


The medical face mask is designed to protect the patient from foreign bodies, covers the mouth and nose, and provides a barrier to minimize the direct transmission of infectious droplets between staff and patients. The mask can also be worn in non-medical areas to reduce the risk of spreading infection. The product is intended for single use only.

## WARNINGS AND LIMITATIONS

- Ensure that you put on the product correctly.
- Dispose the mask at the end of a shift. Do not touch the mask with your hands.
- Wash or disinfect your hands after taking off the mask.
- Wearing the mask does not make you immune to infectious germs.

## APPROVALS AND STANDARDS

- Medical filtering half-mask according to EN 14683, Type IIR
- Bacterial filtration efficiency  $\geq 98\%$

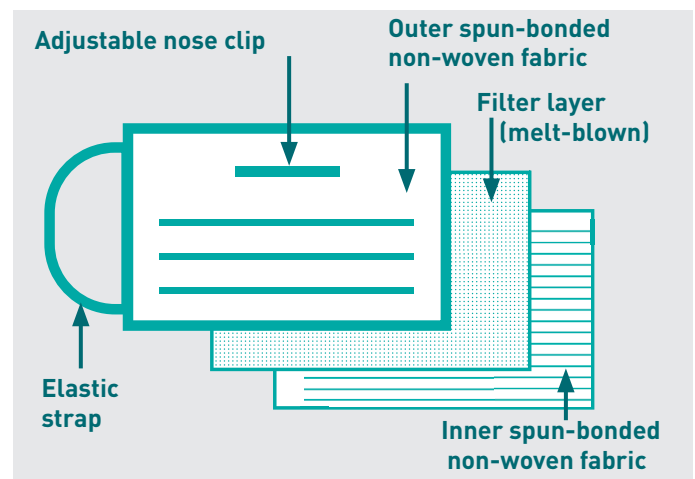

## TIME LIMIT OF USE

Single use, max. one shift

## SHELF LIFE AND STORAGE

- The face mask can be stored for two years in original packaging.
- Ensure that the product is still within its storage period before use. The product should be stored in a clean and dry place within the temperature range of  $-20\text{ }^{\circ}\text{C}$  and  $+40\text{ }^{\circ}\text{C}$  at a maximum relative humidity of 60 %.
- Use the original packaging for storage and transport

